CeraLux

Information and directions for use
Informations concernant l'utilisation

Foglio illustrativo
Gebrauchsinformation

HyCeram[®]
Finest Color Components ■ ■

LuxArt[®]
Unique Color Accents ■ ■

Luxon[®]
Easy Color Shaping ■ ■

invicon
chemical solutions

CeraLux - Table of contents

- 1. Technical description** 5
 - 1.1. Scope of delivery 5
 - 1.2. Device description 5
 - 1.3. Technical data 6
- 2. Safety information** 6
- 3. Setup** 7
- 4. Operation** 8
 - 4.1. Before using for the first time 8
 - 4.2. Light curing 8
 - 4.3. Exposure times 9
 - 4.4. Setting the exposure time 9
 - 4.5. Changing the exposure time 10
 - 4.6. Light output sub-menu – acoustic signal tone on/off 10
- 5. Nitrogen kit** 11
 - 5.1. Description of function 11
 - 5.2. Scope of delivery 11
 - 5.3. Nitrogen quality 11
 - 5.4. Nitrogen kit installation 12
 - 5.5. Nitrogen light curing 16
- 6. Troubleshooting** 17
- 7. Servicing** 17
- 8. Warranty** 18
- 9. EU declaration of conformity** 18

Symbols used

The symbols used in the text are defined below.

- Symbol for a hazard
- Immediate danger to life
- Danger of serious injury
- Danger of minor injury

- Symbol for a hazard
- Risk of material damage
- Risk of damage to the environment

Symbol for useful notes and information

Symbol for required action

Proper use

CeraLux is an industrial light curing unit for curing LuxArt® and Luxon® colour design materials in larger quantities.

It is also suitable for light curing Bond II or Bond THX NF composite media from the HyCeram® colour design system.

1. Technical description

1.1 Product contents

- CeraLux light-curing unit
- Power cable
- Insert plate
- User Information

*optional: object slide tree

1.2. Device description

CeraLux is an industrial light curing unit for curing LuxArt® and Luxon® colour design materials. It is also suitable for light curing Bond II or Bond THX NF composite media from the HyCeram® colour design system. The CeraLux unit cures in the spectral area of 465-475nm with a visible blue light. We recommend optionally the slide tree, maximum 4 pieces find place.

Note!

The unit is not suitable for light-curing materials outside of this spectral area; particularly not for UV-curing materials.

*Slide tree, optionally available

1.3 Technical data

Unit dimensions	300mm/350mm x 377mm x 310mm (Height of front/back x width x depth)
Weight	26kg
Light curing area volume	300mm x 300mm x 170mm
Electrical output	150 watt
Number of LEDs	81 units
LED wavelength area	465 - 475nm
Temperature range	10°C - 30°C
LED dimming range	1 - 99%
Minimum curing time	1 second
Maximum curing time	99 minutes / 59 seconds

2. Safety information

- ▶ Each unit leaves the factory after a final inspection and a test run. However, if the unit is damaged or changed due to transport or other circumstances, it should not be put into operation. Check the unit for damage at regular intervals.
- ▶ All electrical and other supply connections may only be connected in accordance with the technical data. The values must not be exceeded or fallen short of. Observe the VDE regulations. Only use a permitted grounding pin.
- ▶ The regulations regarding any required tests (e.g. in accordance with trade associations or electrical systems) are a matter solely for the buyer/user of the product.
- ▶ This unit is designed for use in accordance with the technical description and therefore must only be used according to its purpose. Any other use is not permitted for safety reasons. In such cases, the manufacturer liability will cease to exist.
- ▶ The unit may only be used by persons who have been instructed in handling the unit and are expressly assigned the task of its operation.
- ▶ Only set up the unit in a clean, heated room with ventilation.
- ▶ The unit must not be covered during operation and free space around the unit of 1m must be adhered to.
- ▶ The surface upon which the unit is set up must be even, straight and stable enough for use.
- ▶ No objects may be placed on top of the unit.

- ▶ If the expected working results are not fulfilled, do not continue to use the unit. The owner is responsible for any consequential damage resulting from the continued use of the unit.
- ▶ To avoid water, fire and explosion damage, block all supplying energy sources after completing the work or after an interruption of work, and secure these against unauthorised access.
- ▶ During servicing, cleaning and repair work as well as before opening the unit, always switch off the main switch and remove the network plug.
- ▶ Repairs may only be carried out by qualified and certified experts.
The system may only be operated with the door closed!

3. Setup

The unit is supplied in a ready-to-use condition and can be operated with a standard 230 V/16 A connection. Observe the safety information for setting up the unit in chapter 2.

Place the unit in your desired working location, easily accessible from all sides.
Make sure you have a clean working environment with as little dust as possible.

4. Operation

4.1 Before using for the first time

Plug the network cable into the connection socket on the back of the unit and connect it with electrical power.

4.2 Light curing

Carry out the light curing of HyCeram Bond II or HyCeram Bond THX NF and LuxArt® /Luxon® colour materials with the following working steps:

Step 1

Open the front door and equip the transparent slide-in surface with the coated objects or optionally the slide tree.

Warning! Cure three-dimensional objects briefly beforehand if necessary with the TwinLux hand-held light curing unit.

Step 2

Insert plate: Place the fitted insert plate on the guide rails provided and close the front door.

Slides: Mount and place the optional slides (max. 4 slides) and close the front door.

Step 3

Set the desired exposure time and activate the exposure by pushing the START button. (for details see point 4.4)

Step 4

At the end of the exposure time, you will hear an acoustic signal tone.

Open the front door and remove the light-cured objects for further processing.

Step 5

A thin, inhibited layer (smear layer) forms on the surface of the colour element during light curing. Roll or wipe off this layer and then further process the objects.

Note:

This smear layer is an effect related to the technology during light curing made from oxygen from the air and has no influence on the strength or quality of the colour elements created.

Note:

If you are using the optionally available nitrogen box, this inhibited layer can be avoided and a hard, glossy surface develops during the exposure process. (see point 5)

4.3 Exposure times

Observe the following exposure times for a sufficient exposure:

Product	Curing time
HyCeram® Bond II	10 minutes
HyCeram® Bond THX NF	30 minutes
LuxArt®	15 minutes
Luxon®	10 minutes

4.4 Setting the exposure time

The desired exposure time can be set as follows:

Step 1

Switch on the main switch on the back of the unit

The screen on the top right of the device displays minutes and seconds.

Step 2

Press the setting button and set the desired exposure minutes by turning the dial. Confirm your selection by pressing the settings button lightly

Step 3

Press the START button to begin the exposure process

4.5 Time correction

You can change the desired exposure time while the unit is in operation as follows:

Step 1

Hold down the START button for 3 seconds. **The time display returns to the original setting**

Step 2

Briefly press the setting dial and make the correct time setting in accordance with point 4.4.

4.6 Light output settings and acoustic signal tone on/off

The CeraLux unit is preset with 100% of the maximum possible light output and with an acoustic signal tone at the end of the exposure time.

If desired, both parameters can be modified as follows:

Step 1

Hold down the settings button for 3 seconds. **The sub-menu for setting the light output appears on the screen.**

Step 2

Set the desired light output by turning the settings button (1-99%) and confirm this by pressing the settings button briefly. **B on-off appears on the screen.**

Step 3

Define on or off by turning and confirm by pressing the settings button briefly. **END appears on the screen.**

Step 4

Confirm the END by pressing the settings button briefly. **The minutes / seconds display appears on the screen.**
The exposure can be continued following point 4.4.

Note:

Invicon recommends the preset light output of 100% for curing times as described in point 4.3.

Note:

The set light output of 100% is the recommended light intensity and should not be changed without a serious reason.
The exposure times stated in point 4.3 refer to the preset light output of 100%.

Warning!

If you increase the light output, you are significantly reducing the lifespan of the LED lights!

5. The nitrogen kit

5.1 Description of function

CeraLux units can be operated with an optionally available nitrogen kit.

This allows you to light-cure in a protective gas atmosphere.

LuxArt® or Luxon® colour elements cured using nitrogen automatically obtain a hard and glossy surface without a smear layer.

The nitrogen kit can be integrated into the CeraLux light-curing unit without the use of tools, and even subsequently.

5.2 Scope of delivery

The CeraLux nitrogen kit contains the following components:

- ▶ Transparent nitrogen curing box with removable lid and tube coupling and Gas ventilation module; dimensions 31cm x 26cm x 5.5cm
- ▶ Flow meter with supporting element
- ▶ Connection tube transparent 80cm
- ▶ Pneumatic tube black 200cm

5.3 Nitrogen quality

Only nitrogen with a purity of > 99.95% is suitable for achieving a high-gloss effect.

5.4 Installation nitrogen kit

Step 1

Remove the cover plugs on the left exterior of the CeraLux device by pressing them lightly from the interior.

Step 2

Remove the plastic fastening nut from the threads of the supporting element.

Step 3

Place the flow meter with wall support on the left side wall of the CeraLux unit.

Step 4

Place the lightly bent base element of the wall support onto the lower edge of the light curing unit and press this upwards until the fastening threads of the wall support fit into the hole of the CeraLux unit side wall.

Step 5

Screw the fastening nut inside of the CeraLux unit onto the coupling threads.

Step 6

Insert one end of the transparent connecting hose into the gas coupling piece on the left inside of the CeraLux device.

Step 7

Slide the plexiglass plate control into the appropriate insertion guides and place the hardening box on top.

Step 8

Push the other end piece onto the coupling piece of the hardening box.

Step 9

Connect the flow meter to the pressure reducing valve of the gas cylinder with the black pneumatic hose.


Note:

The pressure reduction valve is not included in the delivery.

Step 10

Turn the gas flow on the gas container open and adjust the gas flow with the flow meter on the left outside of the CeraLux unit to 5 litres per minute (LPM).

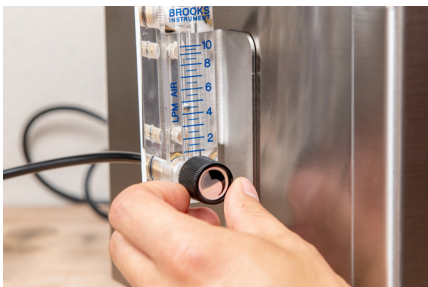Step 11

Set the desired exposure time and activate the exposure by pressing the START button.

Warning!

The gas cannister discharge pressure must not be more than 10 bar.

Note:

To remove the transparent connecting hose again, push the black ring inwards and pull the hose at the same time.

5.5 Nitrogen Light curing

Step 1

Fill the plastic box with the intended objects and close it with the removable lid.

Warning! Allow the pin of the lid to lock into place in the nut of the front wall.

Step 2

Place the filled box onto the guide rails inside the CeraLux.

Warning! Do not place the box on the floor of the CeraLux unit.

Step 3

Start the gas flow at the gas container, adjust the gas throughput with the flow meter on the left exterior of the CeraLux to 8 litres per minute (LPM) and flood the nitrogen box with the gas.

Step 4

Stop the gas flow after 3 minutes. It is recommended to turn off the gas directly at the gas canister to prevent any unnecessary loss of protective gas.

Step 5

Start the exposure. Duration 5 minutes. Unnecessary loss of protective gas.

Note:

If in individual cases, the degree of gloss is not completely achieved after a standard exposure process, increase the exposure time and the gas flow and repeat the exposure process.

6. Troubleshooting

Error	Remedies
The display is off.	Check whether the unit is connected to power. Check whether the unit switch on the back of the unit is at the right position. Check the fuse underneath the unit switch on the back. If the above remedies do not lead to the desired success, contact the manufacturer.
The display is working but the status display is off.	Check whether the time is running on the display.
There is no signal tone at the end of the time.	Check whether the signal tone function has been switched off in the sub-menu. Description according to point 4.6.
The curing degree is not sufficient.	Check whether the time was set according to the specifications under point 4.3. Check the output of the system, description according to point 4.6. Prolong the time until you reach the desired result.

7. Servicing

Clean the insertion plate and nitrogen box carefully on occasion with water and detergent.

The CeraLux is practically maintenance-free. Clean the top of the unit with a dry cloth. The ventilation slits must never be closed!

- 

Warning!

Warning!

Warning!

Do not clean with alcohol!

Never open the unit box! Danger of injury due to electric shock.

Only qualified and certified experts may carry out repairs or replace electrical components.

8. Warranty

The warranty period is 2 years from the date of delivery.

The warranty is only applicable for material or manufacturing faults. The warranty claims do not relate to natural wear, or to damage which occurs after the transfer of risk as a result of incorrect or improper use, excessive use, unsuitable operating materials and such chemical, electrochemical or electrical influences which are not stipulated in the contact.

Improperly carried out changes or repairs and maintenance work carried out by the ordering party or a third party, for example, void any warranty claims.

9. EU Declaration of Conformity

Product name(s)	CeraLux
Manufacturer	Invicon Chemical Solutions GmbH Schweizerstrasse 96 6830 Rankweil AUSTRIA

We confirm under our sole responsibility the conformity of the above mentioned products with the following standards:

Low Voltage Directive 2014/35/EU EN 60598-1:2008 + A11:2009

EMV Directive 2014/30/EU EN 55015:2013
 EN 61547:2009
 EN 61000-3-2 bis -3
 EN 61000-4-2 bis -6, -8, -11

Ecodesign Directive 2009/125/EU VO 1194/2012

RoHS Directive 2011/65/EU

Place of the exhibition	Rankweil
Date of issue	27.05.2019
Name	Peter Bubendorfer
Function	Managing Director

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Invicon chemical solutions GmbH